



Enteral Nutrition Pumps


Model #: _____
Serial #: _____
Date: _____

Provider Name: _____
Phone #: _____
Contact Person: _____

Setting up your equipment:

- Home Use:
 - Fill bag and hang on IV Pole. (Picture 1)
 - Attach tubing to pump. (Picture 2)
 - Pump should be set to OFF before plugging into electrical outlet.
 - Plug in Input 120V, 60Hz, 200 mA
- Ambulatory:
 - Check for battery charge.
 - Fill bag. (Picture 3)
 - Attach to pump.
 - Place into ambulatory bag.

Picture 1

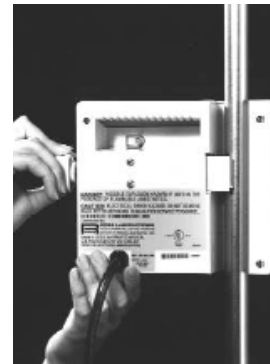

Picture 2



Picture 3



Picture 4




Using your equipment:

- Turn pump dial to the Set Rate position, and set designated food rate (ml/hr). (Picture 4)
 - Locate correct dose and set to instructions
 - Press/Turn/Set to run.
- * If product is not filled at bedside, it is advisable to close slide clamp before transporting.
* Make sure slide clamp is open for set up.



Picture 5

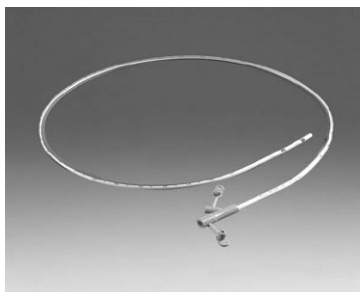

Picture 6



Frequent replacement items:

- Tubing (Picture 5)
- Feeding Bags (Picture 6)

Special Procedures:

- Ambulatory Pump:
 - Must have correct ambulatory bag for pump & tubing bags to work properly. Brands may differ. Check instructions for placement of pump and bag.
 - A new pump should be plugged into an A/C outlet before initial use for directed amount of time to assure full battery charge.

Safety Issues:

- Alarm Sounding:
 - Check for kink(s) in the tubing.
 - Check for empty bag (do not operate pump with empty bag)
 - Check for low battery
 - Check for clogging in the line, slide the clamp to assure it's completely open.
 - Check drip chamber to be sure it is not more than ¼ full or coated with product.
 - Be sure all bags and tubing are compatible with the pump you are using.

Maintenance:

- The outside of the pump can be cleaned with warm soapy water.
- The pump cavity can be cleaned with a cotton swab or soft cloth and warm soapy water. DO NOT CLEAN WITH ALCOHOL.
- Pump should be turned off or recharge battery (if available) when not in use.
- DO NOT SUBMERGE OR AUTOCLAVE the pump or charger.

Please contact your equipment provider for any questions, adjustments, or repairs.

Your local provider:



©Copyright The VGM Group 2002

B.L.U.E. Sheet photos, drawings, and text courtesy of Golden Technologies, Invacare, Leisure-Lift, Permobil, Puritan Bennett, Pride, Respironics, Ross, Sammons Preston, Smith-Nephew, Sunrise Medical